

3 in 1 Pressotherapy Machine



Important Notice

- 1) Please read the user manual before any operation.
- 2) Do not put any liquid onto or into the machine. Be careful about the electric connections.
- 3) Avoid operations on pregnant women, or those people with epilepsy, pacemakers or any serious psychological or medical condition, or infections in the area treated
- 4) Please disassemble or assemble the equipment by professional personnel.
Please check if there is any connection inside the machine dropped when some function does not work.
- 5) Clean and sterilize the spare parts after each operation.

1, Principle and Effects

This aircraft is a most popular multi-function weight-loss device, combined with functions of far-infrared fat dissolving, air-pressure lymphatic drainage, as well as low-frequency stimulation.

1) Far Infrared

The far-infrared slimming treatment can help increase basic metabolic rate (BMR), which is the fastest and natural way to lose weight.

As the body warms, in addition to sweating, the blood thins and the speed at which it circulates around the body increases, while the body tries to cool itself down. The effect is to raise the temperature at the extremities of the body to the same temperature at the core of the body – a temperature at which body fat can be broken down.

2) Pressotherapy


The pressotherapy, with a set of cloth-like accessory and 10 pieces of air bags, which allows for pressure by cycles or all-together for both the whole body and certain parts of the body, is a compression system designed to increase the venous and lymphatic flow and enhance extra-cellular fluid clearance.

3) EMS

By creating electronic pulse to stimulate special points on the body, EMS (Electro Muscle Stimulation) can effectively regulate bio-electric and endocrine system and accelerate metabolism, reaching remarkable effect of weight losing and skin tightening.

2, Panel and Installation Indication

1) Panel

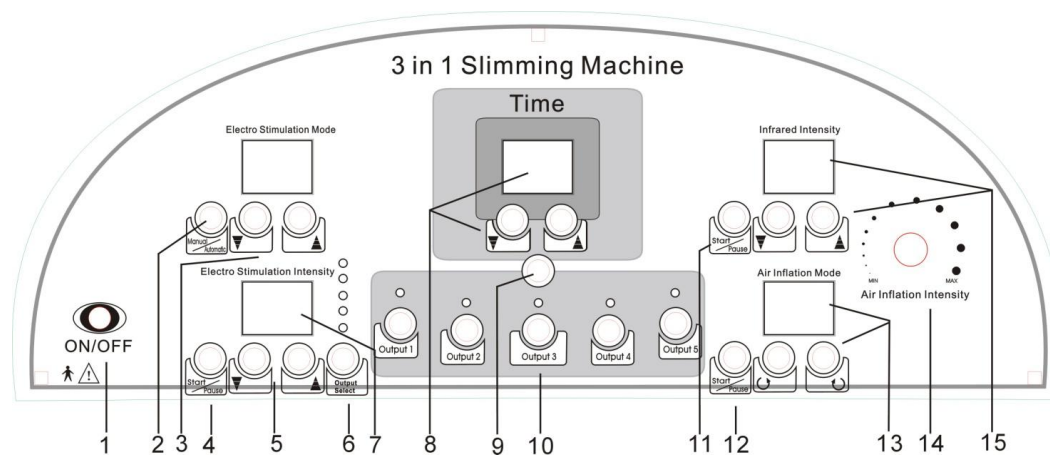

1. Power switch
2. Manual/auto select
3. Wave mode select
4. Start/pause for EMS
5. EMS intensity adjust
6. EMS outlet select
7. EMS intensity display
8. Working time display and adjust
9. Start/pause for selected outlet
10. Outlet select for air inflation
11. Start/pause for infrared
12. Start/pause for air inflation
13. Air inflation mode and adjust
14. Air inflation intensity adjust
15. Infrared intensity and adjust

2) Connector orders on the suit



3) Connectors in the machine

Number for air connector




4) Connections




① EMS cables to be connected to machine

② Far infrared cables and pressotherapy tubes to be connected to machine

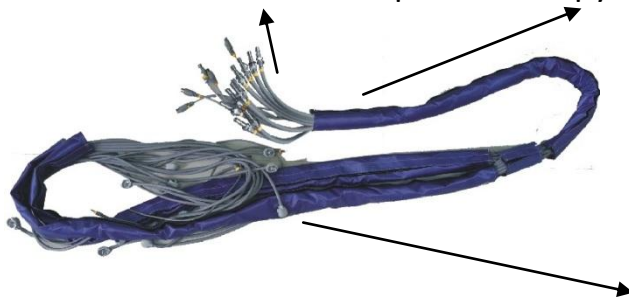

③ Far infrared cables and pressotherapy tubes to be connected to the suit

④ Connection details

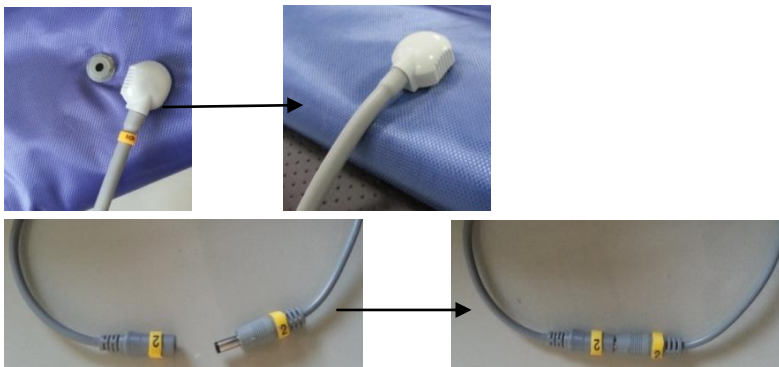

3, Recommended Operations

1) Connect the output port to the slimming clothes and electrode patch respectively.

Here are 2 methods to wear the slimming clothes and electrode patches,

One is to put the electrode patches to the parts with which you are not satisfied at, and put on the slimming clothes for waist, arms and legs. The second is to wrap the preservative film and then put on the slimming clothes, then after finishing infrared treatment take off the preservative film and put the electrode patch on the parts.

2) Pressotherapy

①Select suitable time, and pressotherapy intensity (it is recommended to close the knob from the beginning of operation, then adjust it if necessary, after the cloth is inflated well).

②Select the air inflation mode.

Press the left button to select mode, press the right button to select inflate time.

②-1, There are four modes for air pressure (refer to the number at the left)

Mode 1, it inflates from 0 to 9, it deflates 0 when it inflates 1, and deflates 1 when it inflates 2,,,, then it inflates 0 to 9 together, and the n it repeats the program.

Mode 2, it inflates from 0 to 9, it deflates 0 when it inflates 2, and deflates 1 when it inflates 3,,,, then it inflates 0 to 9 together, and the n it repeats the program.

Mode 3, it inflates from 0 to 9, it keeps inflating 0 when it inflates 1, and keeps inflating 0 and 1 when it inflates 2, and keeps inflating 0, 1 and 2 when it inflates 3, and keeps inflating 0, 1, 2 and 3 when it inflates 4,,,,, then it repeats the program.

Mode 4, the machine inflates at Three lines at the same time,

"0" → "1" → "2" → "3" → "4" → "5"

"6" → "7" → "6" → "7" → "6" → "7"

"9" → "8" → "9" → "8" → "9" → "8"

Then it repeats.

②-2, The numbers at the right means the interval time, 1 to 9 means 1 to 9 seconds.

3) Far infrared operation

- ①Please wrap wet towel before put on the slimming clothes to avoid burning.
- ②Select a proper intensity from low to high according to client's feeling. For each of the 5 outlets, if it is working there will be a little red spot sparkling in the right corner of the display. It will not work if its corresponding connector in the cloth is not connected to the machine.

4) Electro Stimulation

- ①Press button "6" to select outlet, the top one is outlet "5" corresponding to "5" in the connectors. Adjust intensity and press button "4" to start this outlet.

②Wave mode select

Press button “2” to select manual or automatic adjust (in automatic state, the red light will be twinkling)

There are Five modes to choose,

1> Program one: Low Frequency to effect on fat cell and make the fat reduce from the fat cell and then achieve the result of losing weight .Be suitable for any body part.

2> Program Two: Movement for deeply losing weight. It makes the fat move from the deep layer. Be suitable for the deep fat layer

3> Program three: Exercise muscle to losing weight. Be suitable for the beehive-like fat and with the excellent effect of tighten.

4> Program four: Continuous Wave to improve the lymph circulation and reduce the fat from the fat cell .Be suitable for block of lymph circulation with the obvious effect of de-toxin.

5> Program five: Impulse Wave to release the muscle and improve the fat movement.

4, Maintenance for the cloth accessory

1) Do not leave the cloth being heated up without patient inside.

2) Please open the cloth to abstract heat and clean just after finishing the operation.

3) Take care of the connections, do not pull on the pipes roughly, do not bend the pipes excessively.

5, Packing List

Item Name	Unit	Quantity	Remark
Mainframe	set	1	
Cloth accessory and mat	set	1	
Connector	set	1	
Stimulation pad and connector	set	5	
Fuse and power wire	set	1	
Instruction	set	1	

6, Specifications

Voltage	☐ 100-110V ☐ 220-240V		
Frequency	☐ 50HZ ☐ 60HZ	Power	300 W
Package	106*58*48 (machine)	Gross Weight	40 KG